

Pobrano ze strony

<https://desco.pl/>

Chcesz wiedzieć więcej? Skontaktuj się z nami:

☎ [+48 693 919 400](tel:+48693919400)

✉ dh@desco.pl



BEVPOR PW Bottled Water

Filter Cartridges



BEVPOR PW filters ensure the microbiological safety of bottled water whilst protecting the purity and essential characteristics of the source water.

The inert and highly asymmetric PES membrane provides validated microbial retention to industry regulated contaminating organisms. Combined with hydrophilic properties for easy integrity testing, BEVPOR PW filters provide assured performance throughout their service life.

The incorporation of an integral prefilter layer allows graded retention throughout the depth of the filter to resist blockage, resulting in increased capacity and long service lifetime.

BEVPOR PW filters have been designed to provide a cost-effective solution to the microbial stabilization of bottled water by providing increased process control with increased operational efficiency.

Features

Validated retention to industry regulated micro-organisms

Inert materials of construction

Easily integrity tested in-situ

Integral depth prefiltration layer

Benefits


Ensures the safety of the water prior to bottling

Protects the purity and essential characteristics of the source water


Assured filtration performance

Increased throughput to blockage

Performance Characteristics

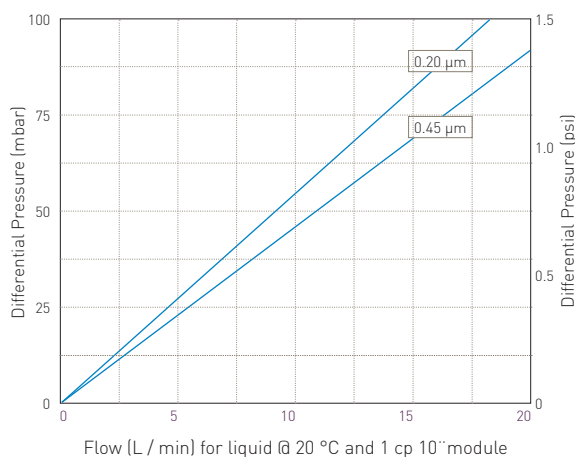

Filtration Stage

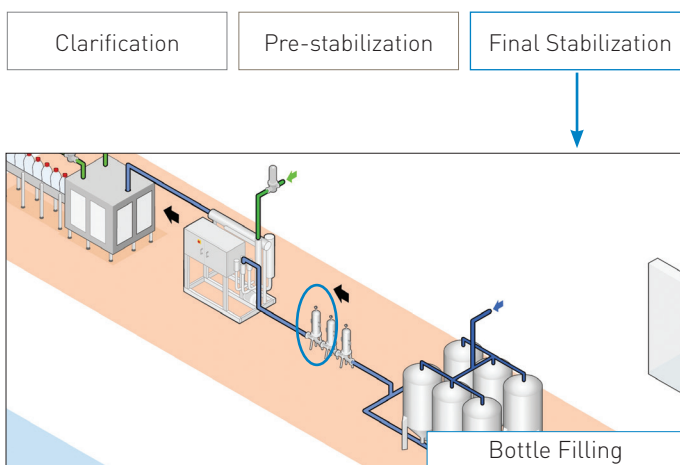

Specifications

Materials of Construction

■ Filtration Membrane:	Polyethersulphone
■ Prefilter Layer:	Polyester
■ Upstream Support:	Polyester
■ Downstream Support:	Polyester
■ Inner Support Core:	Polypropylene
■ Outer Protection Cage:	Polypropylene
■ End Caps:	Nylon
■ End Cap Insert:	316L Stainless Steel
■ O-rings:	Silicone / EPDM

Food Contact Compliance

Materials conform to the relevant requirements of FDA 21 CFR Part 177, current EC1935 / 2004 and current USP Plastics Class VI - 121 °C.



Recommended Operating Conditions

Up to 70 °C (158 °F) continuous operating temperature and higher short-term temperatures during CIP to the following limits:

Temperature		Max Forward dP	
°C	°F	(bar)	(psil)
20	68	5.0	72.5
40	104	4.0	58.0
60	140	3.0	43.5
80	176	2.0	29.0
90	194	1.0	14.5
>100 (steam)	>212 (steam)	0.3	4.0

Effective Filtration Area (EFA)

10" (250 mm) Up to 0.6 m² (6.45 ft²)

Cleaning and Sterilization

BEVPOR PW cartridges can be repeatedly steam sterilized in-situ or autoclaved at up to 130 °C (266 °F). They can be sanitized with hot water at up to 90 °C (194 °F) and are compatible with a wide range of chemicals. Please refer to our Clean-in-Place support guide or contact your local Parker representative for more information.

Retention Characteristics

0.2µm BEVPOR PW filters have been validated to provide sterile effluent after bacterial challenge testing following ASTM F838-05 methodology on 10" cartridges with more than 10⁷cfu per 10" cartridge using *Brevundimonas diminuta*.

In addition, challenges with the following EU regulated organisms have been performed.

Organism	LRV when challenged with a minimum of 10 ⁷ cfu per cm ²	
	0.20	0.45
<i>Serratia marcescens</i>	FR	FR
<i>Escherichia coli</i>	FR	FR
<i>Enterococcus faecalis</i>	FR	FR
<i>Clostridium perfringens</i>	FR	FR
<i>Pseudomonas aeruginosa</i>	FR	9.1
<i>Brevundimonas diminuta</i>	5	-

*FR - Fully retentive during challenge

When expressed as titre reduction "FR" equates to >10⁷ per 10" module.

Integrity Test Data

All filters are flushed with pharmaceutical grade purified water prior to despatch. They are integrity tested to the following limits:

Diffusional Flow Test Parameters	Micron Rating	
	0.20	0.45
Test Pressure (barg)	1.7	1.4
Test Pressure (psig)	25.0	20.0
Max Diffusional Flow per 10" (ml /min)	16.0	16.0

Manufacturing Traceability

Each filter cartridge displays the product name, product code and lot number. Additionally, each module displays a unique serial number providing full manufacturing traceability.

Ordering information

BPW		-			-		A	
Code Length (Nominal)		Code Micron		Code End Cap (10 inch)		Code O-rings		
1	10" (250 mm)	02	0.20 µm	C	Fin / 226 Bayonet	S*	Silicone	
2	20" (500 mm)	04	0.45 µm	D	Fin / 222	E	EPDM	
3	30" (750 mm)			G	Recess / 222			
4	40" (1000 mm)			R	BF / 222 Bayonet			

*Silicone O-rings supplied as standard

*Silicone O-rings supplied as standard

PHP & PHL
HOUSING RANGE
AVAILABLE